

# nEXT TURBOMOLECULAR PUMPS



nEXT  
240



## INNOVATION AND RELIABILITY

Edwards nEXT is the ultimate experience in turbomolecular pumps. nEXT turbomolecular pumps are built on decades of experience and are drawing from our tried and trusted EXT and STP ranges. nEXT pumps offer superior performance, reliability and end user serviceability, setting the benchmark for scientific turbomolecular pumps.

## **ADVANCED TECHNOLOGY**

nEXT has been designed to combine all the latest technological advances in turbomolecular pumps with some new thinking in design for manufacture, delivering a truly class leading product.

The nEXT platform brings a high level of modularity to offer maximum flexibility for customer application and requirements. Each pump is available in different internal configurations to offer differing functionality and performance.

Our nEXT pumps come in different variants "D" Duplex comes with both turbomolecular and drag stages for improved tolerance to higher backing line pressures. The "D" variants offer superior pumping speed and compression across all gas species. Triplex "T" variants feature turbomolecular, drag and Edwards unique regenerative pumping stages for the ultimate in compression ratio and boost technology for unique vacuum system rationalisation.

The "H" variant has been physically tuned to offer an improvement where an application has focus on light gas compression.

### **Superior performance**

Exceptional pumping speeds and compression ratios

## **Proven reliability for peace of mind**

Huge install base of turbo pumps

### **Flexible solutions**

Bespoke design service available

### **Ease of use**

Integrated intelligent controls

### **Extended lifetime and low cost of ownership**

Fully end user serviceable

### **Large variety of standard variants**

Enhanced customer choice

# APPLICATIONS & FLANGE SIZES

## Applications


	nEXT55	nEXT85	nEXT240	nEXT300	nEXT400	nEXT730	nEXT930	nEXT1230
<b>Research &amp; Development</b>								
Chamber evacuation, coating systems, turbomolecular pump systems	✓	✓	✓	✓	✓	✓	✓	✓
<b>High Energy Physics</b>								
Beam lines, accelerators, mobile pump carts, turbomolecular pump backing, laser evacuation, medical systems	✓	✓	✓	✓	✓	✓	✓	✓
<b>Mass Spectrometry</b>								
GCMS, LCMS, ICPMS, MALDI, inorganic MS, RGA, surface science, leak detectors	✓	✓	✓	✓	✓	✓		
<b>Electron Microscopy</b>								
TEM, SEM, EPMA, SPM sample prep benches	✓	✓	✓	✓	✓			
<b>Industrial</b>								
Glove boxes, coating systems, XRD/XRF systems, leak testing, energy, furnaces, medical technologies	✓	✓	✓	✓	✓	✓	✓	✓

You can be assured Edwards has the application expertise and vacuum solution to meet your needs.

## Inlet flange sizes

	NW40	DN63 ISO-K	DN63-ISO CF	DN100 ISO-K	DN100 ISO CF	DN160 ISO-K	DN160 ISO CF	DN200 ISO-K	DN200 ISO CF	DN200 ISO-F
nEXT55D	x*	x	x	x*						
nEXT85D	nEXT85H	x	x	x						
nEXT240D	nEXT240T			x	x					
nEXT300D	nEXT300T			x	x					
nEXT400D	nEXT400T					x	x			
nEXT730D						x	x	x	x	
nEXT930D								x	x	
nEXT1230H								x	x	X

\*Available by special order

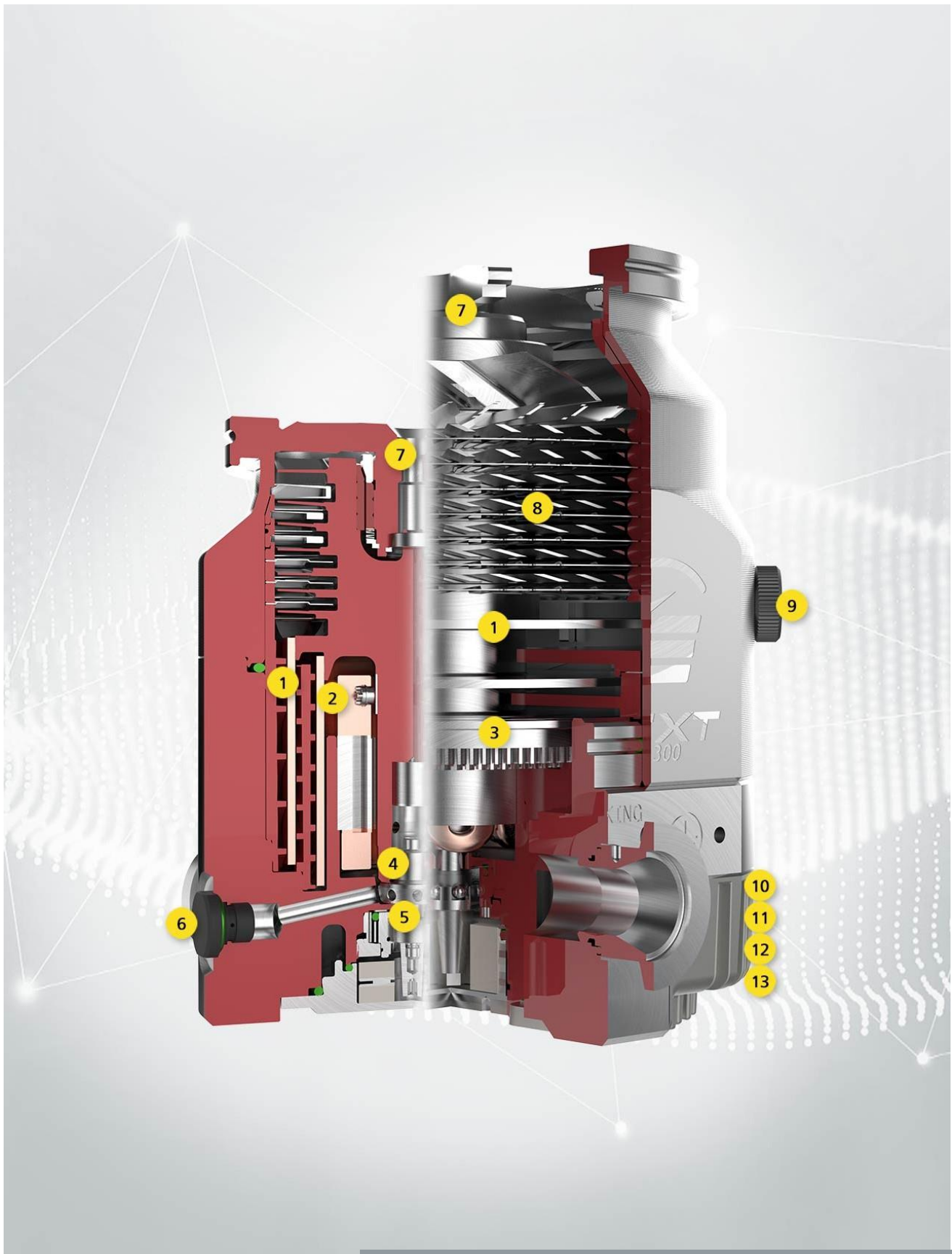

PERFORMANCE YOU CAN  
RELY ON

## **1 - Multiple drag stages**

For higher compression ratios and greater backing pressure tolerance

## **2 - Direct temperature measurements within the pump**

Enables wider operating envelope (nEXT55-nEXT85)

## **3 - Third regenerative stage**

Enables boost port options and higher compression (nEXT240, 300, 400)

## **4 - Patented bearing suspension system**

Minimises noise and vibration transmitted to vacuum (nEXT240, 300, 400)

## **5 - User replaceable oil cartridge**

For a speedy interim service and bearing cartridge for a quick low cost scheduled overhaul (nEXT55-400)

## **6 - Purge port**

For safe operation of pump with specialised gases

## **7 - Permanent magnet upper bearing**

For a hydrocarbon free vacuum, reduced vibration and minimum wear

## **8 - Advanced simulation tool**

Derived rotor design to give better speed and compression performance

## **9 - Manual vent port**

With automatic valve accessories for rapid venting and quick cycle times

## **10 - Manual as well as serial setting of standby speed**

For flexible system tuning

## **11 - Simple parallel operation**

Or more sophisticated serial control in both RS232 and RS485 protocols

## **12 - Automatic wide operating voltage range**

From 24V to 48V d.c. for versatility in system integration (nEXT730-1230 48V)

## **13 - Sensorless drive**

For high efficiency and compactness with reduced stray magnetic fields





## nEXT55 AND nEXT85

The most compact pumps of the nEXT range with a significantly reduced height and improved performance in a smaller package.

These pumps offer pumping speeds of 55 l/s and 85 l/s for nitrogen, they provide a high pumping density, greater than other pumps in its class, with almost double the pumping speed of similar sized turbo pumps.

**The rotor has been designed to optimise pump performance and achieve both higher speeds and higher compression while maintaining high levels of reliability and low risk to adoptors.**


As with all the pumps in the nEXT range field maintenance is only required every 4 years of operation, and this including replacement of the bearing can be performed by the end user.

nEXT55 and nEXT85 bring with them the benefits in flexibility with comprehensive communication and control options available, as well as a full set of accessories, the ideal choice in deployable instruments or portable applications where a compact footprint or lower weight are key factors.




# DIMENSIONS & PERFORMANCE


nEXT55

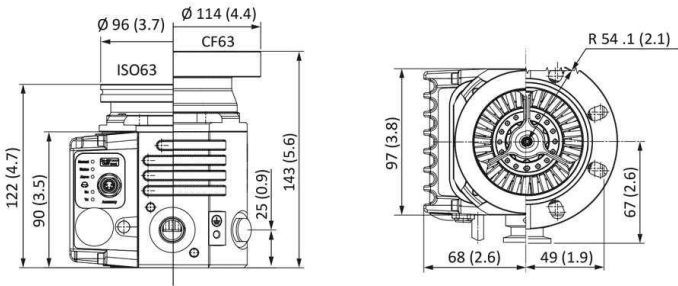


nEXT85

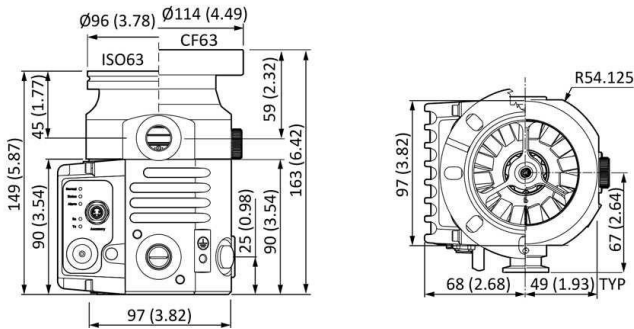


nEXT55D

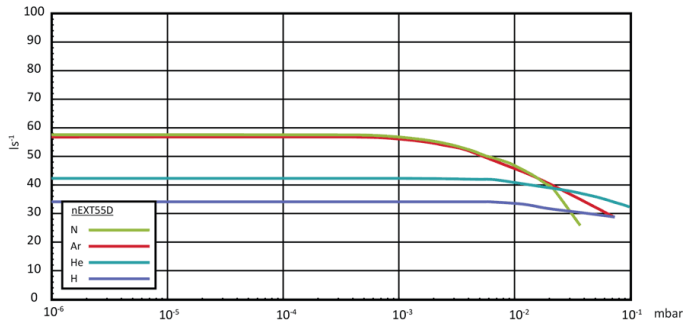

nEXT85 DN63

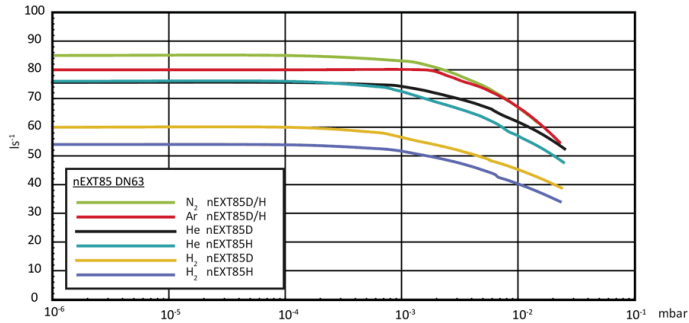

nEXT85 DN40

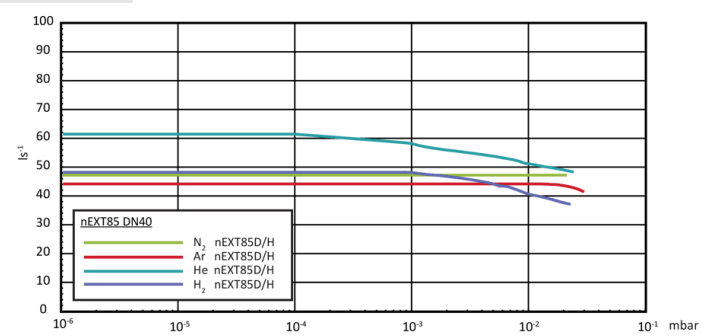



## nEXT240, nEXT300 & nEXT400

The innovative pumps, designed to provide high vacuum performance in a compact size.

## **A wide variety of sizes**

Giving OEMs and end users a greater choice and the flexibility to tailor the most appropriate vacuum solution to meet their individual needs.

## **Each is available in two variants**

- the D-Type combines turbo and drag stages;
- the T-Type adds Edwards unique fluid dynamic stages and the option of additional booster ports for increased system rationalisation.

See Boost technology custom interface splitflow variants are also available in three sizes for further vacuum system optimisation.

## **Field replaceable oil cartridge and bearing assembly**

The pumps feature a field replaceable oil cartridge and bearing assembly and the user is notified as to when service intervention is required.

This enables easy maintenance as users can carry out servicing in-house, which reduces the cost of ownership.

## **Low power consumption and standby mode**

The efficient pumps have low power consumption and also feature a standby mode, which allow users to make further energy savings.

nEXT240-300-400 pumps have extra low vibration and low magnetic field signature variants for sensitive electron microscope applications.

## **Compact design**

The compact design means the pumps fit together neatly in high product density applications.

They are easy to configure and have an intelligent control interface accepting a 24 to 48Vdc input power supply voltage range and can be controlled using a simple parallel control or serial communications in both RS232 or RS485 formats.

# BOOST TECHNOLOGY

For our Analytical Instrument Customers the “T” variants of nEXT240,300,400 employ a unique regenerative pumping stage which can be used to ‘boost’ the performance of the backing pump.

## BOOST TECHNOLOGY


**A much simplified vacuum solution with greatly improved pumping speeds and system power reduction**

The examples below shows how boost can be used to either reduce the size of the primary pump or increase the gas flow into the instrument on a differentially pumped system.


It also illustrates how the number of turbo pumps required can be reduced from two discrete pumps to a single splitflow pump with two inlets to achieve the same vacuum performance.

Customers in general laboratory and R&D applications will also benefit from the improved compression achieved with the "T" variant.

To take maximum advantage of boost technology, please contact Edwards.


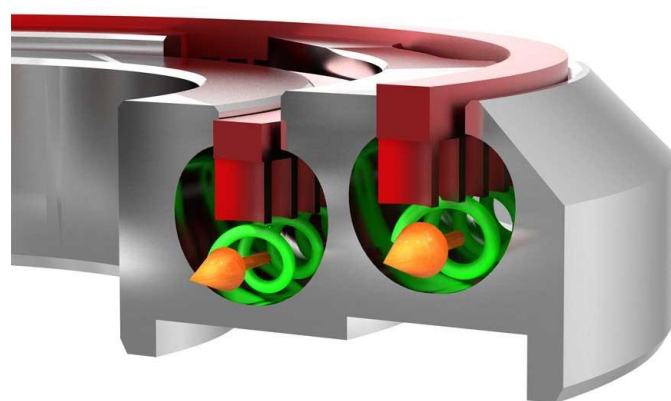


### Regenerative Stage

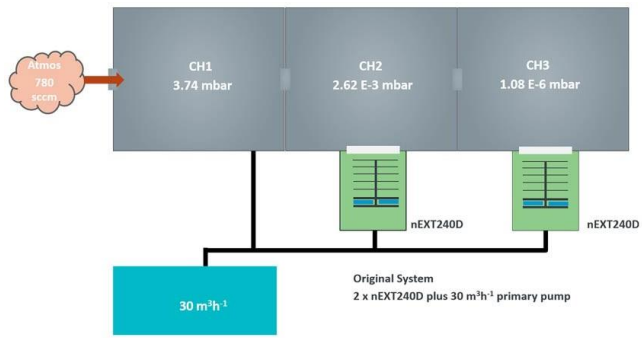

## Original System


## Boost Option 1

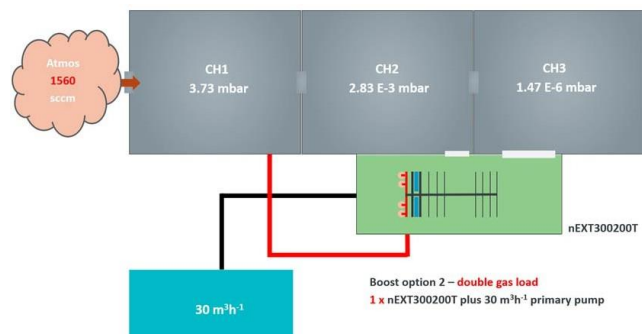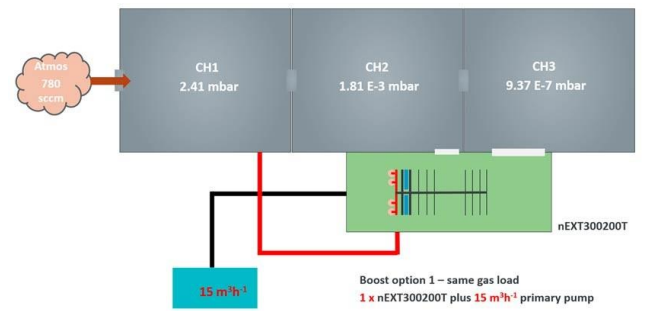

## Boost Option 2




## Screenshot taken from TransCalc HSM

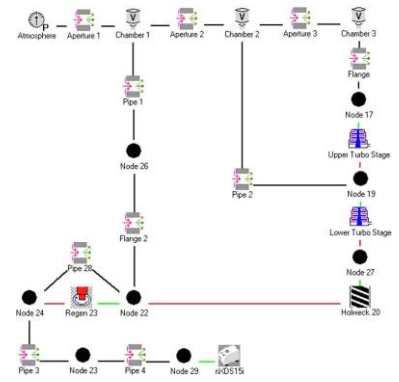


# DIMENSIONS & PERFORMANCE


nEXT240D

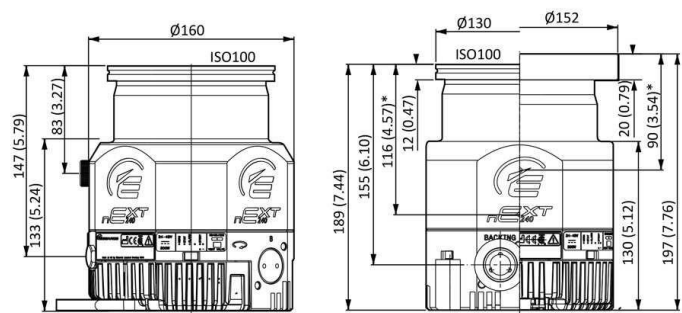


nEXT240D

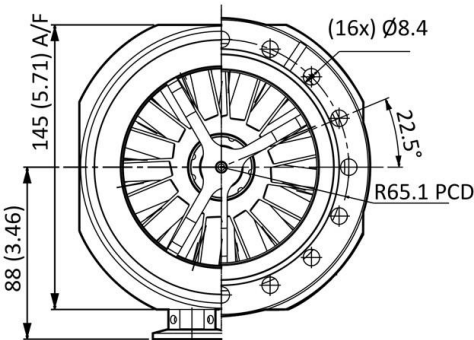


nEXT300D

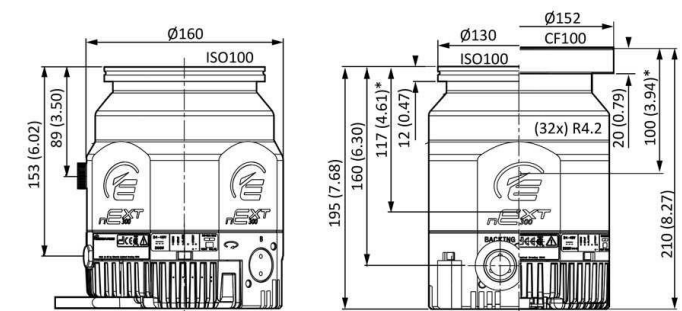


nEXT300D

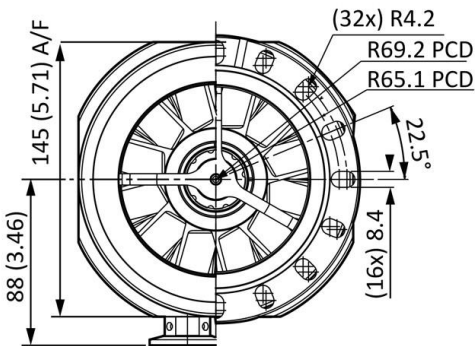


nEXT400D

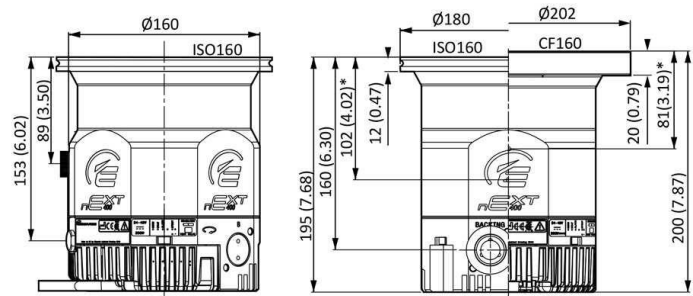


nEXT400D

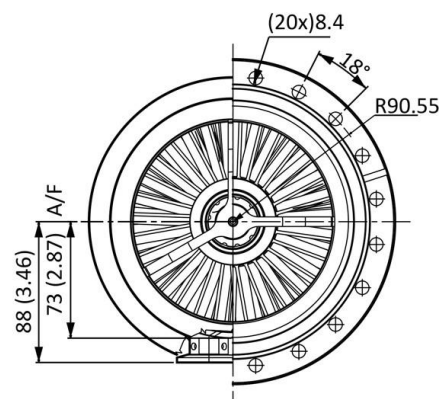


nEXT240

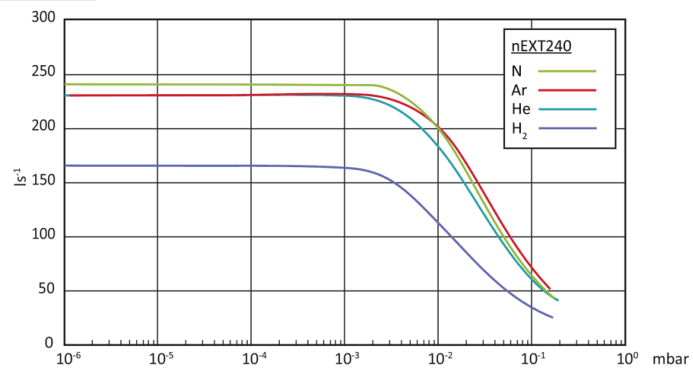

nEXT300

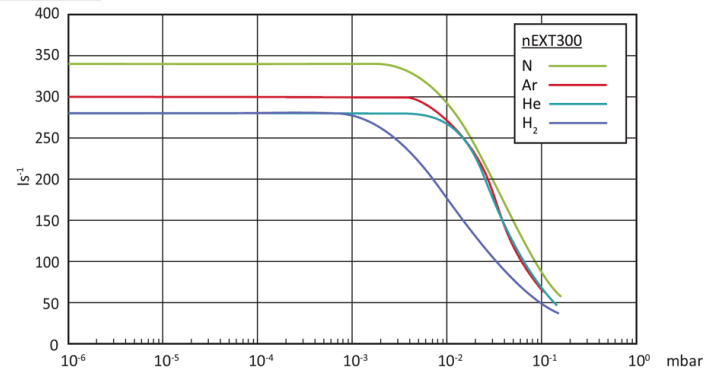

nEXT400

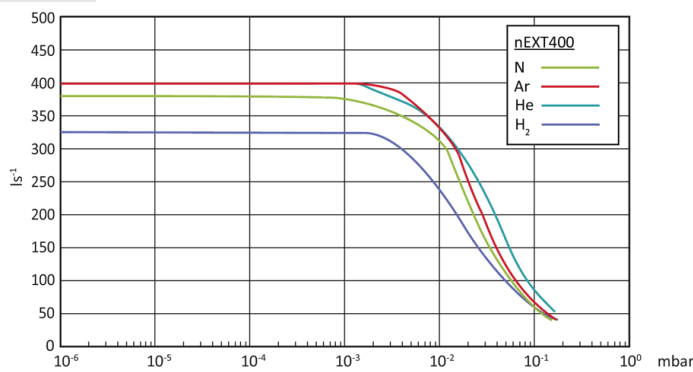



## nEXT730, nEXT930 & nEXT1230

Larger pumps offering nitrogen pumping speeds of 730, 925 and 1250 l/s respectively.

## **Designed to meet a wide range of requirements**

As well as addressing the general R&D market, where high pumping speeds are often required, these pumps are also designed to meet the requirements of the coating market and other diffuse market sectors such as:

- Heat treatment
- Furnace applications
- E-beam welding
- Etch
- Ion implant
- Degassing
- Cylinder evacuation

## **We match the specific applications of OEM**

For our OEM customers, derivative versions of these products including split flow variants can be developed to match specific applications.

## **Operate in any orientation**


These compact pumps are able to operate in any orientation, and are supported by a full range of accessories for cooling, venting, powering and control.

## **The pumps feature bearings with a typical life time of at least 4 years with no maintenance.**


When required, a bearing replacement can be undertaken by the customer themselves or they can take advantage of our other service support options.

# DIMENSIONS & PERFORMANCE


nEXT730D

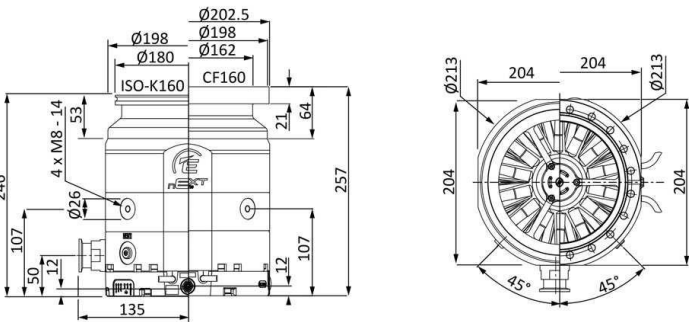


nEXT930D

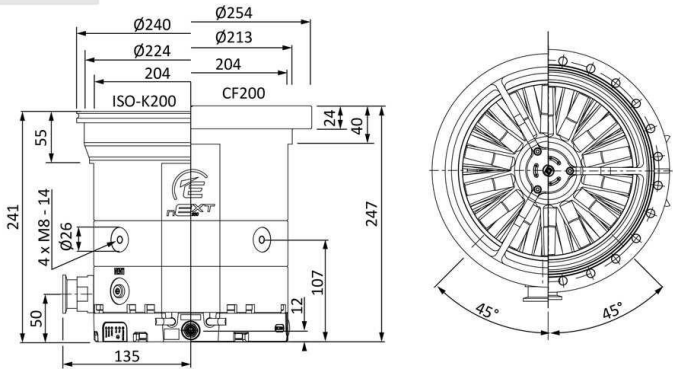

nEXT1230H

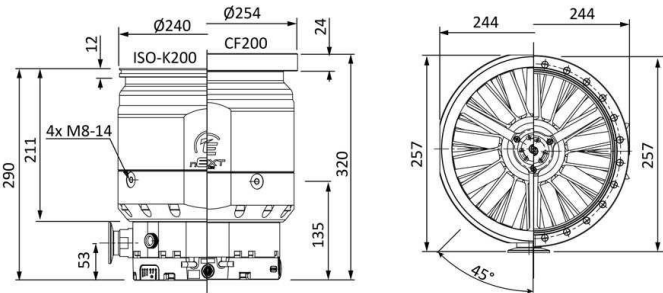


nEXT730D

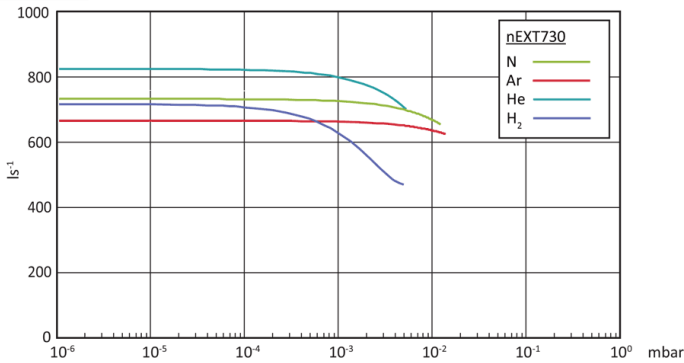




nEXT930D

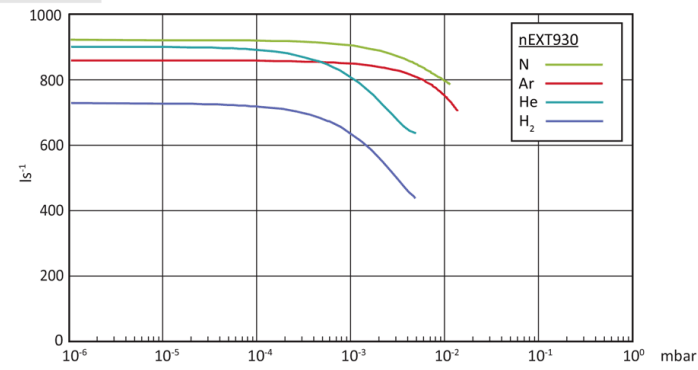

nEXT1230H

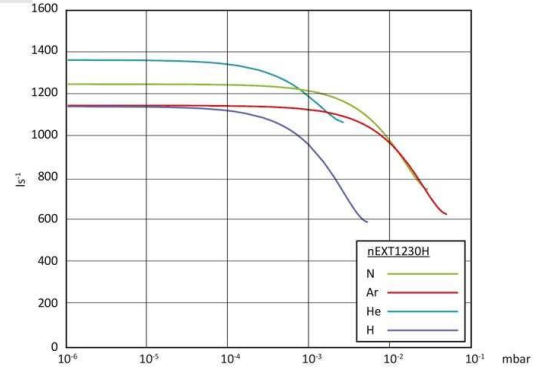

# CREATING THE PERFECT SOLUTION

With Edwards, you can quickly and accurately create the perfect vacuum solution for your needs.



## A partnership approach to system design

**When an off the shelf pump will not meet your requirements for space or performance, our Bespoke Product Development (BPD) team will develop a customised vacuum solution to turn your requirement into reality.**

At Edwards, it all starts with a vacuum expert gathering your application details!

We want to understand what pressures you wish to achieve, what gas flows you have, how much space is available etc.

Our expert then uses a number of in-house and publicly available modelling tools at his disposal to optimise your vacuum system. Developed by our Analytical Services group, these tools are used to model complete vacuum systems from atmosphere down to ultra-high vacuum (UHV). This software has been developed to give rapid simulation of the behaviour of the proposed vacuum solution to ensure that it perfectly meets your requirements.

Accurate computer modelling offers you the chance to streamline your development cycle, avoiding a costly iterative approach and delivering a quicker time to market. Please contact Edwards to take advantage of this service.

## TIC

### Turbo and Instrument Controller

Automatically recognises and supports one turbomolecular pump from the nEXT range, one backing pump (nXDS or diaphragm) plus three Edwards active gauges on top of that cooling and vent valve support is provided directly from the controller. Backing pump power is provided for a compact 24V diaphragm pump (on 200W versions only), or where greater pumping speeds are required, nXDS and large XDS pumps can be controlled directly via the backing pump connector on a TIC 200 mains backing pumps (up to RV12) may be controlled via an optional relay box.

The relay box can also be used to control a mains heater band and backing line isolation valve. Time delays and normal speed signals may be used to control events such as turbo start and there is a comprehensive selection of protection and safety interlock features. The TIC turbo controller may be either rack or bench mounted and provides a useful hub for the flexible operation of a wide range of vacuum system configurations.



## TAG

### Turbo and Active Gauge

A small, compact, low cost pumping system controller, which is suitable for a wide range of vacuum applications. It is a 24V controller that is compatible with all Edwards nEXT turbomolecular pumps. In addition to a turbomolecular pump it can control a backing pump, a vent valve, an air cooler and an Edwards active gauge. The TAG is controlled by an easy to use interface. A large clear LED display shows the pump speed or vacuum pressure. The compact size of the controller is ideal for use on benchtops or suitable mobile platforms.



# SERVICE

## End user serviceability

**New technologies employed in nEXT have enabled the pumps to be serviced by the end user in the field.**


An oil lubrication cartridge change can be performed where fitted, typically in less than 5 minutes. A full bearing change can also be performed by the end user in around 10 minutes on all nEXT pumps. Both with the minimum of specialist tooling. These simple interventions will, in many cases, mean that the pump never requires a full return to base service during its lifetime.

nEXT turbomolecular pumps will advise the user when a service is due and what level of intervention is required.

The user is alerted to a service request by a simple flashing LED sequence on the pumps and by serial comms notification.

Flexibility is again key as these simple services can be performed either by the end user, on site by an Edwards field service technician, or the pump can be returned to an Edwards service hub.

Using remote diagnostics, a user can interrogate the pump to determine how long it is to the next service so that a proactive approach to preventative maintenance can be planned.

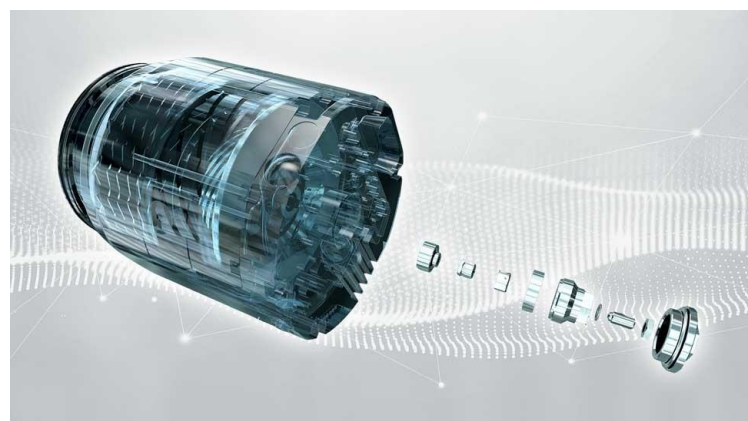



## — Prolonged peace of mind

**Extending the new equipment warranty gives you a simple opportunity to add peace of mind to your purchase of new equipment**

Should a fault occur as a result of a manufacturing defect, equipment is expressly repaired or replaced.

Cover is available on many of our products allowing the original factory warranty to be extended from 12 months to 2 years and beyond.

## — Your global partner

**We understand the importance of local support.**

Edwards has a number of major service facilities located throughout the world, each location is supported by an extensive team of engineers and technicians to provide local, rapid response and great value service.

All our service operations are conducted at the highest international standards in accordance with ISO9001 (Quality), ISO14001 (Environmental), and OHSAS18001 (Workplace safety).



# TECHNICAL SPECIFICATIONS

## nEXT55, nEXT85

		nEXT55	nEXT85		
Inlet flange		DN63 ISO-K or DN63 CF	NW40	DN63 ISO-K or DN63 CF	DN100 ISO-K
Inlet pumping speed ls <sup>-1</sup>	N <sub>2</sub>	55	47	84	86
	Ar	55	44	80	84
	He	41	61	78	80/78 (D/H)
	H2	27	49/44 (D/H)	60/54 (D/H)	60/54 (D/H)
Compression ratio (D)	N <sub>2</sub> /Ar	>1 x 10 <sup>11</sup>			
	He	6.9 x 10 <sup>5</sup>	8 x 10 <sup>6</sup>		
	H <sub>2</sub>	2.9 x 10 <sup>4</sup>	2 x 10 <sup>5</sup>		
Compression ratio (T/H)	N <sub>2</sub> /Ar	NA	>1 x 10 <sup>11</sup>		
	He	NA	2 x 10 <sup>7</sup>		
	H <sub>2</sub>	NA	5 x 10 <sup>5</sup>		
Backing/interstage/boost ports		NW16			
Vent/purge port		1/8" BSPP			
Critical backing pressure (D/H)	mbar	18			
Critical backing pressure (T)	mbar	NA			
Bake out water cooled/forced air cooled max.	°C	120/115°			
Recommended backing pump*		nXDS6i			
Normal rotational speed (rpm)		90,000			
Start time to 90% speed (sec) D/H (T)		90			
Mass (kg) D/H (T)	ISO	2.47	2.9	3.0	3.2
	CF	3.5		4.4	

\*A smaller backing pump may be used depending on application.



## nEXT240, nEXT300, nEXT400

		nEXT240	nEXT300	nEXT400
Inlet flange		DN100 ISO-K or DN100 CF	DN100 ISO-K or DN100 CF	DN160 ISO-K or DN160 CF
Inlet pumping speed ls <sup>-1</sup>	N <sub>2</sub>	240	300	400
	Ar	230	280	380
	He	230	340	390
	H <sub>2</sub>	165	280	325
Compression ratio (D)	N <sub>2</sub> /Ar	>1 x 10 <sup>11</sup>	>1 x 10 <sup>11</sup>	>1 x 10 <sup>11</sup>
	He	3 x 10 <sup>5</sup>	1 x 10 <sup>6</sup>	1 x 10 <sup>8</sup>
	H <sub>2</sub>	1 x 10 <sup>4</sup>	5 x 10 <sup>4</sup>	5 x 10 <sup>5</sup>
Compression ratio (T/H)	N <sub>2</sub> /Ar	>1 x 10 <sup>11</sup>	>1 x 10 <sup>11</sup>	>1 x 10 <sup>11</sup>
	He	1 x 10 <sup>6</sup>	3 x 10 <sup>6</sup>	>1 x 10 <sup>8</sup>
	H <sub>2</sub>	1.5 x 10 <sup>4</sup>	1 x 10 <sup>5</sup>	1 x 10 <sup>6</sup>
Backing/interstage/boost ports		NW25	NW25	NW25
Vent/purge port		1/8" BSPP	1/8" BSPP	1/8" BSPP
Critical backing pressure (D/H)	mbar	9.5	9.5	10
Critical backing pressure (T)	mbar	20	20	20
Bake out water cooled/forced air cooled max.	°C	120/115°	120/115°	120/115°
Recommended backing pump*		RV12/nXDS10i	RV12/nXDS10i	RV12/nXDS10i
Normal rotational speed (rpm)		60,000	60,000	60,000
Start time to 90% speed (sec) D/H (T)		115 (150)	145 (190)	180 (210)
Mass (kg) D/H (T)	ISO	5.7 (6)	5.7 (6)	6.5 (6.8)
	CF	8.8 (9.1)	8.5 (8.8)	9.5 (9.8)

\* A smaller backing pump may be used depending on application.

## nEXT730, nEXT930, nEXT1230

		nEXT730D		nEXT930D		nEXT1230H		
Inlet flange		DN 160 ISO-K	DN 160 CF	DN 200 ISO-K	DN 200 CF	DN 200 CF	DN 200 ISO-F	DN 200 ISO-K
Inlet pumping speed ls <sup>-1</sup>	N <sub>2</sub>	730		925			1250	
	Ar	665		865			1150	
	He	820		905			1350	
	H <sub>2</sub>	715		735			1150	
Gas throughput mbar ls <sup>-1</sup>	N <sub>2</sub>	14		14			12	

## nEXT730, nEXT930, nEXT1230

		nEXT730D		nEXT930D		nEXT1230H		
Gas throughput mbar ls <sup>-1</sup>	Ar	3.5		3.5		4		
	He	21		21		>20		
	H <sub>2</sub>	>>14		>>14		>20		
Compression ratio	N <sub>2</sub>	>1 x 10 <sup>11</sup>						
	Ar	>1 x 10 <sup>11</sup>						
	He	1.2 x 10 <sup>8</sup>		1.2 x 10 <sup>8</sup>		4 x 10 <sup>+8</sup>		
	H <sub>2</sub>	4 x 10 <sup>6</sup>		4 x 10 <sup>6</sup>		1 x 10 <sup>+7</sup>		
Ultimate pressure with 2-stage oil sealed rotary vane pump ISO-K/CF	mbar	<3.5 x 10 <sup>-9</sup>	<6 x 10 <sup>-10</sup>	<3.5 x 10 <sup>-9</sup>	<6 x 10 <sup>-10</sup>	<5 x 10 <sup>-10</sup>	<5 x 10 <sup>-9</sup>	
Backing/interstage/boost ports		NW40						
Normal rotational speed	rpm	49 200		49 200		42 000		
Start time to 90% speed (sec) D/H (T)	min	2.5		2.5		3		
Cooling water consumption	l/h	60						
Critical backing pressure	mbar	15						
Mass (kg) D/H (T)	kg	14.6	19.6	15.4	21.7	32.6	24.9	23.7
Recommended backing pump*		nXRi, XDS35i, E2M28*						
Bake out water cooled/forced air cooled max.	°C	n/a	100	n/a	100	100	n/a	n/a
Vent/purge port		G1/8"						

\*Please contact your local representative to discuss the correct option for your application.

# ORDERING INFORMATION

## Pumps

Product description	Order number
nEXT55D NW40 NW16 80W	B8E210A01
nEXT55D CF63 NW16 80W	B8E210C01
nEXT55D ISO63 NW16 80W	B8E210B01
nEXT55D ISO100 NW16 80W	B8E210101
nEXT85D ISO63 NW16 80W	B8G210B01
nEXT85D CF63 NW16 80W	B8G210C01
nEXT85D ISO100 NW16 80W	B8G210101
nEXT85D NW40 NW16 80W	B8G210A01
nEXT85D ISO100 NW25 80W	B8G240101
nEXT85iD ISO63 NW16/16 80W	B8G211B01
nEXT85iD ISO63 NW16/25 80W	B8G214B01
nEXT85H ISO63 NW16 80W	B8G410B01
nEXT85H CF63 NW16 80W	B8G410C01
nEXT85H NW40 NW16 80W	B8G410A01
nEXT85iH CF63 NW16/16 80W	B8G411C01
nEXT85H ISO100 NW16 80W	B8G410101
nEXT240D ISO-K100 160W	B81200100
nEXT240D CF100 160W	B81200200
nEXT240T ISO-K100 160W	B81300100
nEXT240T CF100 160W	B81300200
nEXT300D ISO-K100 160W	B82200100
nEXT300D CF100 160W	B82200200
nEXT300T ISO-K100 160W	B82300100
nEXT300T CF100 160W	B82300200
nEXT400D ISO-K160 160W	B83200300
nEXT400D CF160 160W	B83200400
nEXT400T ISO-K160 160W	B83300300

## Pumps

Product description	Order number
nEXT400T CF160 160W	B83300400
nEXT730D ISO-K160 NW25	B8J200300
nEXT730D CF160 NW25	B8J200400
nEXT930D ISO-K200 NW25	B8K200D00
nEXT930D CF200 NW25	B8K200F00
nEXT1230H CF200 NW40	B8N4A0F00
nEXT1230H ISO-F200 NW40	B8N4A0E00
nEXT1230H ISO-K200 NW40	B8N4A0D00
nEXT1230H CF200 NW40 INV	B8N4A0FU0
nEXT1230H ISO-F200 NW40 INV	B8N4A0EU0
nEXT1230H ISO-K200 NW40 INV	B8N4A0DU0

Other interstage port positions available upon request

## Accessories and spares nEXT55/85

Pump	Product description	Order number
Controller <sup>(1)</sup>	TAG controller	D39592000
	TAG power supply	D39592800
	TIC 200 turbo and instrument controller	D39721000
Cooling	WCX85 water cooling kit (4 position)	B8G200833
	ACX85 air cooler connector fitted	B8G200820
Venting	N/O TAV5 vent valve connector fitted	B8G200834
	N/C TAV5 vent valve connector fitted	B8G200835
Bakeout	CF63 flange heater 110 V	B8G200823
	CF63 flange heater 240 V	B8G200824
Service	Oil cartridge kit	B8G200828
	Bearing and oil cartridge kit	B8G200811
	Bearing replacement tool kit	B8G200845
Miscellaneous	Accessory Y adaptor	B8G200837
	Accessory cable 90 degree/extension	B8G200836
	Accessory connector bare wired	B8G200839
	nEXT85/EXT75DX base mounting adaptor	B8G200838

(1) Denotes need second annotation nEXT730 and bigger pumps need their own power supply required accessory. Others optional depending on application.

## Accessories and spares nEXT240/300/400

Pump	Product description	Order number
Controller <sup>(1)</sup>	TAG controller	D39592800
	TAG power supply	D39592800
	TIC100 turbo and instrument controller	D39722000
Cooling	nEXT radial air cooler	B58053175
	nEXT axial air cooler	B58053185
	nEXT water cooler	B80000815
Bakeout	CF100 flange heater 100-120 V	B58052773
	CF100 flange heater 200-240 V	B58052774
	CF160 flange heater 100-120 V	B58052775
	CF160 flange heater 200-240 V	B58052776
Venting	TAV5 solenoid operated vent valve	B58066010
Service	Oil cartridge tool kit	B80000812
	Bearing tool kit	B80000805
	Oil cartridge	B80000811
	Bearing and oil cartridge	B80000810

(1) Denotes need second annotation nEXT730 and bigger pumps need their own power supply required accessory. Others optional depending on application.

## Accessories and spares nEXT730/930/1230

Pump	Product description	Order number
Controllers	TAG controller	D39592000
	TAG power supply	D39592800
	TIC100 turbo and instrument controller	D39721000
Cooling	Air cooling radial nEXT730/930	B8J200800
	Air cooling radial nEXT1230	B8J200801
	Water cooling nEXT730/930, 1/4 inch	B8J200820
Venting	N/O TAV5 vent valve connector fitted	B8G200834
	N/C TAV5 vent valve connector fitted	B8G200835
	Vent port adaptor	B58066011
Inlet screens	Center ring w. prot. screen DN200 ISO-K coarse	B8J200807
	Center ring w. prot. screen DN200 ISO-K fine	B8J200808
	Coarse inlet screen DN 200 CF	B8J200809
	Fine inlet screen DN 200 CF	B8J200810

## Accessories and spares nEXT730/930/1230

Pump	Product description	Order number
Inlet screens	CF160 coarse inlet screen	B80000823
	CF160 fine inlet screen	B80000824
	ISO160 coarse inlet screen	B80000825
	ISO160 fine inlet screen	B80000826
Service	Bearing replacement kit	B8J200827
	Bearing replacement tooling	B8J200845
Mains input cable	Mains input cable 3m EU	B8J200812
	Mains input cable 3m US	B8J200813
	Mains input cable 3m UK	B8J200814
Power supply to pump	nEXT 3m open end cable	B8J200816
	nEXT 5m open end cable	B8J200817
	EPS 800	B8J200819
	nEXT 3m extension cable for EPS 800	B8J200824
	nEXT 5m extension cable for the EPS 800	B8J200825
	2.5m EU EPS 800, IP54 protected	B8J200829
	2.5m US NEMA 6-15P 250V EPS 800, IP54 protected	B8J200830
	2.5m UK EPS 800, IP54 protected	B8J200831
	EPS 800 mount kit to place the power supply directly on the pump	B8J200832
Miscellaneous	Accessory cable 90 degree/extension	B8G200836
	Accessory Y adaptor	B8G200837
	Accessory connector bare wired	B8G200839

(1) Denotes need second annotation nEXT730 and bigger pumps need their own power supply required accessory. Others optional depending on application.

## Accessories and spares

Pump	Product description	Order number
Extension cables	1 m pump to controller cable	D39700835
	3 m pump to controller cable	D39700836
	5 m pump to controller cable	D39700837
Power cables <sup>(-1)</sup>	2 m electrical supply cable UK plug	D40013025
	2 m electrical supply cable EU plug	D40013030
	2 m electrical supply cable US plug	D40013120
Miscellaneous	Vent port adaptor	B58066011



## Accessories and spares

Pump	Product description	Order number
Miscellaneous	PRX10 purge restrictor	B58065001
	TIC relay	D39711805

(1) Denotes need second annotation nEXT730 and bigger pumps need their own power supply required accessory. Others optional depending on application.

Official Distributor in Australia



**EZZI VISION**

*Vacuum and Thin Film Technology*

### CONTACT US

T: 1800 GO EZZI

E: [sales@ezzivision.com.au](mailto:sales@ezzivision.com.au)

W: [ezzivision.com.au](http://ezzivision.com.au)

**VIC:** 13/62 Ramset Drive, Chirnside Park,  
VIC 3116, Australia

**NSW:** Unit 1, 80 O'Riordan St, Alexandria,  
NSW 2015, Australia

**WA:** Unit 11, 24 Baile Road, Canning Vale,  
WA 6155 Australia



3601 0091 01 - November 2021. All rights reserved.  
Edwards and the Edwards logo are trademarks of  
Edwards Limited. Whilst we make every effort to ensure  
that we accurately describe our products and services,  
we give no guarantee as to the accuracy or  
completeness of any information. Edwards Ltd,  
registered in England and Wales No. 6124750,  
registered office: Innovation Drive, Burgess Hill, West  
Sussex, RH15 9TW, UK.